



DISTEL

— c —

DOSING SYSTEM

A 1 litre bottle of DISTEL DOSING SYSTEM concentrate solution generates 2000 measured doses of surface disinfectant. With the DISTEL DOSING SYSTEM every dose is the same volume as the previous one.

The DISTEL DOSING SYSTEM is intended for the cleaning and disinfection of non-invasive medical devices and hard surfaces within the hospital environment.

The DISTEL DOSING SYSTEM is effective against bacteria, yeast and viruses including **coronaviruses from 3 minutes**.

It is also destroys DNA and RNA.



The DISTEL DOSING SYSTEM is CE marked as a Class IIa Medical Device.
 *Dose onto a new wipe. For full instructions please see the DISTEL DOSING SYSTEM User Guide.

Designed to be repurposed and reused, the DISTEL DOSING SYSTEM bottles keep the environment in mind.

A 1 litre bottle of DISTEL DOSING SYSTEM concentrate delivers 10 x 1 litre bottles of working solution which equates to 200 x 5ml doses per bottle.

The DISTEL DOSING SYSTEM can be used with any wipe to create freshly made surface disinfectant wipes on the go. Choose a sustainably sourced wipe to reduce single-use plastic waste.

ELIMINATE SINGLE-USE PLASTIC WIPES

1 LITRE BOTTLE OF CONCENTRATE

1




BOTTLES OF WORKING SOLUTION

10



FRESH DOSES OF SURFACE DISINFECTANT

2000

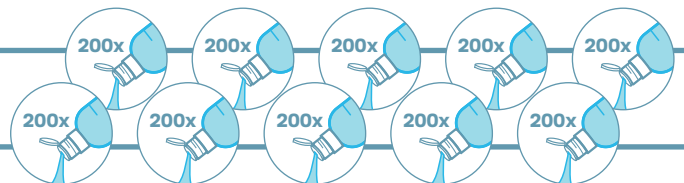

USE THE WIPE OF YOUR CHOICE



MANUFACTURED IN THE UK



REUSABLE BOTTLES



DESTROY CORONAVIRUSES IN 3 MINUTES

THE DISTEL DOSING SYSTEM IS INTENDED FOR THE CLEANING AND DISINFECTION OF NON-INVASIVE MEDICAL DEVICES AND HARD SURFACES WITHIN THE HOSPITAL ENVIRONMENT SUCH AS:




**WORKSTATIONS
(MONITORS
& KEYBOARDS)**



**TROLLEYS
(I.E. SMALL EQUIPMENT
& RESUSCITATION)**



HOISTS

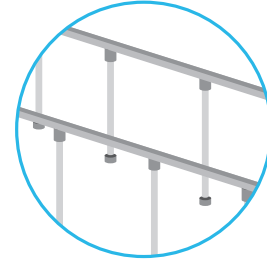

HANDRAILS



A MEASURED 5ML DOSE EVERY TIME

The DISTEL DOSING SYSTEM heads make certain that each fresh dose of surface disinfectant is the same volume. Simply invert the dosing bottle once to prime and a second time to apply 5ml of surface disinfectant.

Each dose is consistent, resulting in uniform disinfection for increased compliance.



**1 LITRE DISTEL DS STARTER PACK
CCH050101**

1 x 1 Litre Bottle of Concentrate - 5 x 1 Litre Dosing Bottles
6 x 5ml Dosing Heads - 1 x Working Solution Bottle Label



**1 LITRE DISTEL DS
CCH050201**

1 x 1 Litre Bottle of Concentrate - 1 x 5ml Dosing
Head - 1 x Working Solution Bottle Label



**1 LITRE DISTEL DS DOSING BOTTLES
CCH050301**

6 x 1 Litre Dosing Bottles - 6 x 5ml Dosing Heads

1 LITRE DISTEL DOSING SYSTEM

**2000 DOSES
PER 1 LITRE
CONCENTRATE
BOTTLE**



**500ml DISTEL DS STARTER PACK
CCH050401**

1 x 500ml Bottle of Concentrate - 5 x 1 Litre Dosing Bottles
6 x 5ml Dosing Heads - 1 x Working Solution Bottle Label



**500ml DISTEL DS
CCH050501**

1 x 500ml Bottle of Concentrate - 1 x 5ml Dosing
Head - 1 x Working Solution Bottle Label



**500ml DISTEL DS DOSING BOTTLES
CCH050601**

6 x 500ml Dosing Bottles - 6 x 5ml Dosing Heads

COMING SOON: 500ml DISTEL DOSING SYSTEM

**1000 DOSES
PER 500ml
CONCENTRATE
BOTTLE**